Universal Chamber Luigs & Neumann



Bath Chamber for Universal Usage


The universal bath chamber allows the use of object holders (26x76), NUNC dishes, Falcon dishes (up to 60mm) and the LN bath chamber inserts. Additionally the chamber can be used with the LN bridge 500, the LN bath chamber manipulator system and the original microscope tables of microscope suppliers. The design of the chamber allows the direct attachement of the reference electrode and in-/out-flow using magnets on the steel surfaces.

Accessories

- Magnets for the attachment of reference electrode, capillaries or perfusion systems
- Bath sensor for measure the temperature of the bath, incl. 1 piece magnet bracket type 2
- o Reference-Electrode, incl. silver chloride wire
- Angle guide-rail for not temp. controlled chambers (for mounting on the LN Mini25/23 Manipulator/Unit)
- Guide-rail for not temp. controlled chambers (for mounting on the LN Mini23 Manipulator/Unit)
- Bath chamber fixture for not temp. controlled chambers (for mounting on the Bridge 500)
- LN inserts for the slice mini and universal bath chamber with groove for PE-Perfusion tubing; black teflon coated, with different drillings (types A to E)
- PE-Perfusion tubing, useable for pre-heating of the ringer medium at the LN insert type D and E





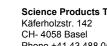

in collaboration with:



Email: biomedical-instruments@science-products.com